


The RADSampler is a compact yet highly sensitive laboratory gamma spectrometer

Developed as a desktop unit, the RADSampler quickly and accurately analyzes any small object to ensure it is free of radioactive contaminants. The RADSampler's lead lined analyzer chamber virtually eliminates ambient background radiation producing ultra-clean, highly accurate test results.

- Compact footprint for desktop application
- Configurable scanning or dwell times
- Available in Basic and Advanced versions
- Easy to operate 5 button keypad
- Remote service and calibration capability

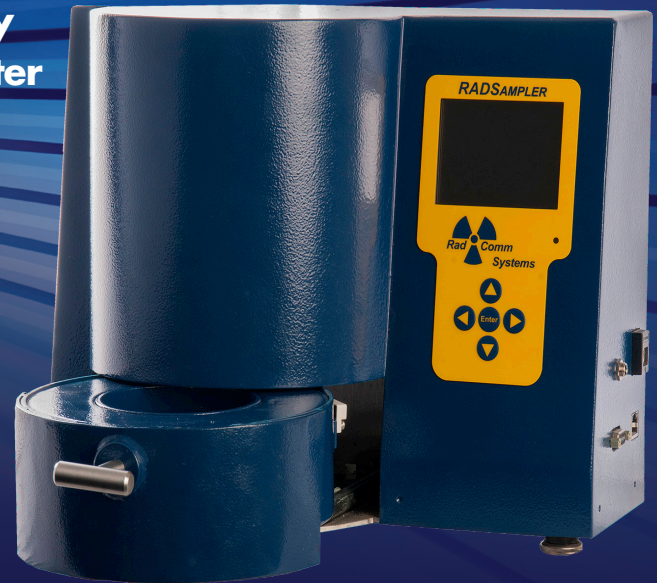

RADSAMPLER

LABORATORY GAMMA SPECTROMETER

Easily configurable to isolate results

RADSampler Basic version allows a user to isolate up to four (4) Regions Of Interest (ROI) of the energy range spectrum and set specific alarm thresholds. This isolation allows for quicker results to detect the four most common isotopes found in the metal manufacturing industries.

RADSampler Advanced takes the user to the next level of detection and identification capability. With a large Isotope Library it can accurately identify specific and multiple isotopes within any given sample. RADSampler Advanced is suitable for any industry that is concerned about potential contaminated product and/or incoming material that requires testing.

Simplified easy to understand results

The RADSampler utilizes powerful Windows based software which allows for uploading of stored spectra results to any PC via USB. The energy results from the ROI windows are automatically summarized in a user friendly report in which the operator can enter information, save, print and/or email the report. Further analysis can be performed by using RadComm's PC Analysis software to view a specific ROI in a detailed histogram format.

Desktop sized laboratory with powerful response

Measuring approximately 8½" (21cm) wide x 12¾" (32cm) long and weighing only 86lbs (39kg), it easily sits on any desktop or industrial shelf. No longer is it necessary to dedicate valuable floor space to accommodate much larger and heavier spectrometers.

RADSampler's large Thallium Doped Sodium Iodide, NaI(Tl), crystal coupled with customer enabled ROI configuration quickly optimizes test results specifically for Am²⁴¹, Cs¹³⁷, Ir¹⁹², and Co⁶⁰. RADSampler Advanced monitors all ROI's providing users with an even greater range of isotopes and analysis.

Detector case consists of:

- Detector dimensions: 9"H (23cm) x 12.69"W (32cm) x 8.5"D (21.6cm)
- Detector case: Painted Aluminum
- Shielding material: Lead
- Detector weight: 86lbs (39kg)
- Maximum sample well size: 3.25" (8.25cm)

Electronics

- Integrated PMT with electromagnetic shielding
- Audio & visual alarm
- Stable low noise, high voltage power supply
- USB connector for both power and data transfer
- Internal Lithium Ion battery for backup purposes
- Internal memory storing up to 900 spectra
- 5 button actuated keypad

Environmental

- Operating temperature: -20°C (-4°F) to +60°C (+140°F)

Display

- Backlit RGB colour LCD 320x240 resolution
- Viewing area: 3.5" (8.9cm) TFT LCD

Software

- Menu driven user interface
- Windows based PC software for data management
- Remote service and calibration capability
- RadView software for integration


Response and Sensitivity

- Energy range: 30 KeV to 3.0 MeV (Gamma)
- RADSampler Basic optimized for Am²⁴¹, Cs¹³⁷, Ir¹⁹², and Co⁶⁰
- Gamma spectrum: 1024 channel
- Detection material: Thallium Doped Sodium Iodide crystal, NaI(Tl)
- Crystal size: 1.5" (38mm) x 2" (51mm)
- Energy resolution: 8.0% or better for 662 KeV
- Calculates dose rate (RADSampler Advanced only)

RadSampler RS00000 Report

Item	Details	Notes	
Time	11:14 AM	Time from RadSampler (hh/mm/ss)	
Date	05/03/2013	Date from RadSampler (dd/mm/yyyy)	
Location	RadComm	Entered by Operator	
Operator Name	John Smith	Entered by Operator	
Heat ID	ABCD123	Entered by Operator	
Sample ID	123	Entered by Operator	
Acquisition Time	10	In Seconds	
R.O.I. #s	CPS	Threshold	Result
R.O.I. # 1	9	13	Passed
R.O.I. # 2	2	5	Passed
R.O.I. # 3	0	5	Passed
R.O.I. # 4	0	5	Passed
Total Counts	49	63	Passed

Sample Result **Passed**

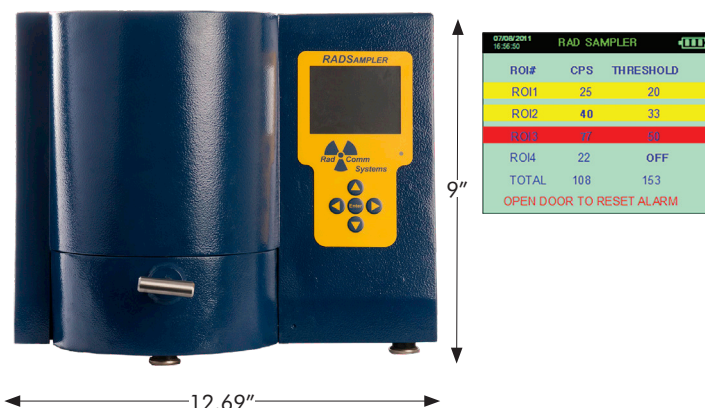

Corporate Head Office

Scott Aikin
Saikin@radcommsystems.com
Jeff Adams
jadams@radcommsystems.com
2931 Portland Drive
Oakville, ON Canada L6H 5S4
Tel. +1 (905) 829-8290
Toll Free. 1 (800) 588-5229
Fax. +1 (905) 829-1406

USA

Joshua Hunter
jhunter@radcommsystems.com
602 E. Lincolnway Ave.
Valparaiso, IN USA 46383
Tel. +1 (773) 680-8430
Toll Free. 1 (800) 588-5229
Fax. +1 (219) 510-5764

Europe

Wim van Hove
wim.van.hove@radcomm europe.com
Watertorenweg 32, 2230 Herselt, Belgium
Tel. +32.14.75.02.13
Fax. +32.14.75.02.16

Black Sea

Osman Tureyyen
metkim@metkim.com
Kavacik, Onarimli Sokak
No 9 Altay Han
Daire: 8 Beykoz, Istanbul, Turkey
Tel. +90 (216) 425-1412
Fax. +90 (216) 425-1432

Latin America

Eduardo Ballesteros
analiticacontrol@prodginy.net.mx
Amatista No.27, El Pedregal
Atizapan de Zaragoza, 52948, Estado De Mexico
Tel. (52) 55 5077 4633/ (52) 55 3615 4910

India

Neelakshi Bhargava
neelakshi@radcommsystems.co.in
C-34, Ground Floor, Sawan Park, Ashok Vihar
New Delhi 110052, India
Tel. +91-9717671924
Fax. +91-9818650118